



Communicator 2.0

Getting Started

An emergency response device with **Intelligent Voice Activation™**



Table Model



Wall Model





Installing the Communicator

1

Open the Communicator device



Hold here and pull off the back cover

2

Press the Pairing button on the Care@Home™ control panel



Control Panel

Pairing Button

3

Insert the batteries and wait for the control panel to beep



4

Close the Communicator



Tamper switch

Tamper pin

WALL MODEL

- 1 Clean and dry the surface
- 2 Peel the strips off the mounting tape and press the Communicator into place

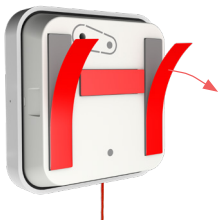


TABLE MODEL

- 1 Place in a convenient place, like a living room table or bedroom night stand.



To mount on a wall, remove the table mount and follow the wall model instructions.




Testing the installation

To test, press the emergency button, or call out the trigger phrase, or pull the cord.

Table Model



Wall Model

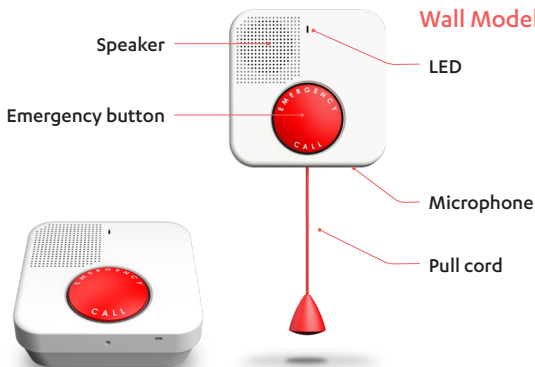

In an emergency



Call out the
trigger phrase



Press
Emergency Call



Pull the cord,
if included

For more information, see: <https://www.essencesmartcare.com/assistance-center>

Legal notice: <https://www.essencesmartcare.com/privacy-notice>

<https://www.essencesmartcare.com> • info@essence-grp.com

© 2023 All rights reserved to Essence SmartCare Ltd.

